


Hidden Analytical Ltd Hidden Analytical Ltd Launches Advanced Software Package for QGA 2.0: MASsoft Professional and QGA 2 Software

Hidden Analytical Ltd is proud to announce the release of its advanced software package, combining MASsoft Professional and QGA 2 software, designed to complement the QGA 2.0 gas analyser. This software suite represents a significant leap in gas analysis technology, offering unparalleled precision and ease of use.

MASsoft Professional: A Multilevel Control Hub

MASsoft Professional is a comprehensive software package allowing for both straightforward control of mass spectrometer parameters and sophisticated manipulation of data. It is integral to all Hidden instruments, including the QGA 2.0, providing a user-friendly interface and efficient data management. Key features include:


- **Ionisation Energy Control:** For soft ionisation of complex mixtures.
- **Data Export:** Compatibility with NIST MS database and external software like Excel, Origin.
- **Device Control:** Manages external devices like MFCs, gas valves, and furnace controllers.
- **Real-Time Data Management:** Allows subtraction of overlapping peaks and outputs data in various formats including percentage and ppm.

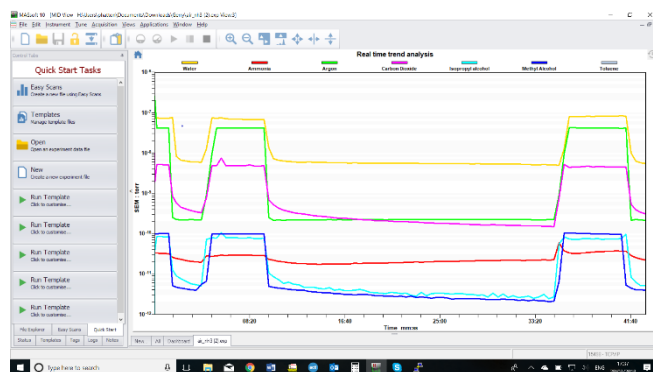

QGA 2 Software: Tailored for Real-Time Gas and Vapour Analysis

The QGA 2 software is a specialized package for quantitative gas and vapour analysis, enabling real-time continuous analysis of up to 32 species with concentrations ranging from 0.1 ppm to 100%. This software offers:

- **Flexible Stream Modes:** Capable of single or multi-stream analysis with up to 80 streams.
- **Intuitive Calibration Routines:** Simplifies both cracking pattern and Relative Sensitivity (RS) measurement.
- **Advanced Analysis Features:** Utilizes template setup routines with automatic spectral removal algorithms for precise quantitative data.
- **Integration with External Devices:** Enhances versatility for a wide range of applications.

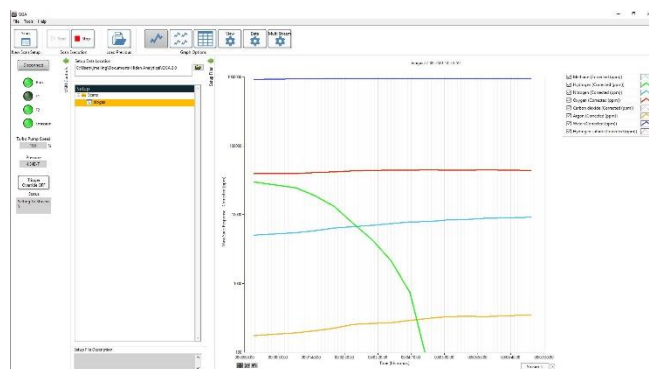

Together, MASsoft Professional and QGA 2 software offer an unmatched level of control, precision, and flexibility, making the QGA 2.0 not just a gas analyser, but a complete analysis solution.

For more information about MASsoft Professional and QGA 2 software, please visit [Hidden Analytical Ltd's website](http://HiddenAnalytical.com) and [QGA 2 software page](#), or send us an email to info@hidden.co.uk.

--- ends ---